
Venen Engel

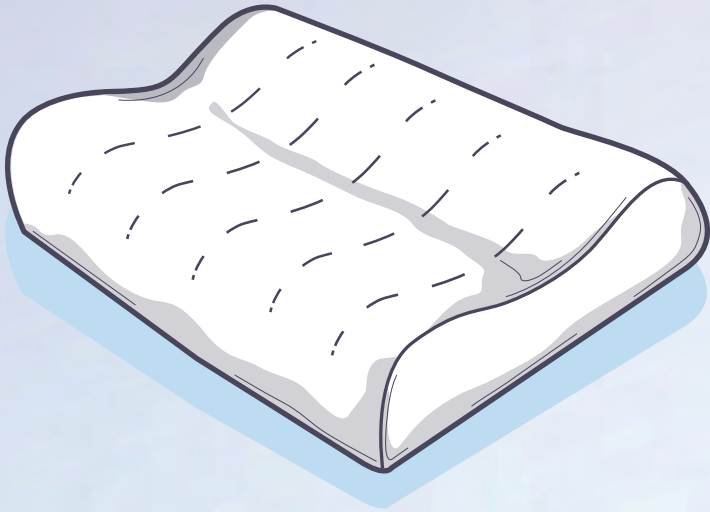VITAL FLOW MASSAGEKISSEN

VITAL FLOW MASSAGE PILLOW

VITAL FLOW MASSAGEKISSEN

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
Web: www.Venen-Engel.de

VITAL FLOW MASSAGE PILLOW


READ ALL INSTRUCTIONS BEFORE USING THE PRODUCT

PACKAGING CONTENT

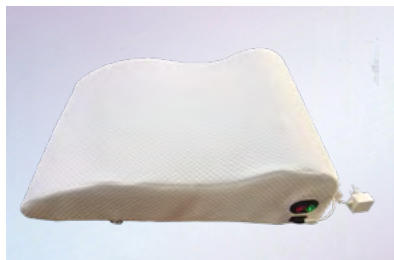Vital Flow Massage Pillow
Power cable

VITAL FLOW MASSAGE PILLOW

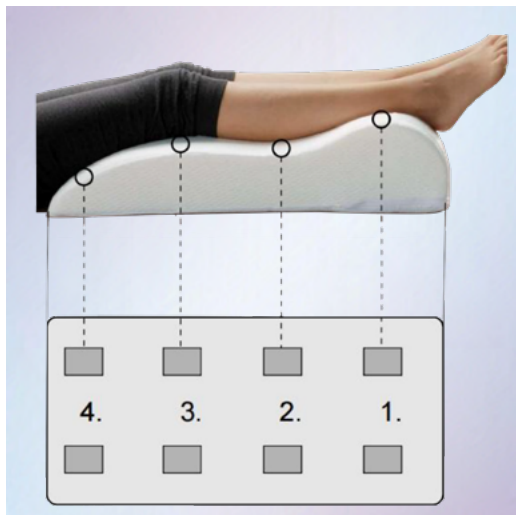Massage points on the pillow

CONTROLLER & CONNECTOR

PRODUCT FEATURES

- 8 massage motors
- 5 massage modes
- 2 heating zones with 3 heating levels
- Intelligent wireless Bluetooth control
- Four massage positions
- Two working modes: power supply via battery or adapter

SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be taken, including the following:

1. This product is for use on a 5V/3A adapter. Make sure that the pillow is connected to an outlet having the same configuration as the plug.
2. Always unplug the pillow from the power source after use and before cleaning.
3. Do not reach for the pillow if it has fallen into water. Disconnect the power cable immediately.
4. An electrical appliance should never be left unattended when plugged in. Unplug from outlet when not in use and/or before attaching or detaching product parts.
5. Do not use outdoors.
6. Do not use the supplied cord as a handle when carrying the product.
7. Keep the cord away from heated surfaces.
8. Close supervision is necessary when the pillow is used by/on/near children or disabled persons.
9. This product is not a toy. Children should not use it unsupervised or use it as a toy.
10. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
11. Do not use near eyes or other highly sensitive areas.
12. To disconnect, push the „**MASSAGE**“  and „**HEATING**“  buttons until the light turns off. Make sure that all functions have stopped, then remove the plug from the outlet.
13. Use the pillow only for its intended purpose as described in this manual. Do not use any attachments that were not originally provided with the unit.

HEALTH INFORMATION

- If you have any concerns regarding your health, consult your doctor before using this product.
- Individuals with pacemakers should not use this product.
- If you experience pain in a muscle or joint for a prolonged period, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.

A heating pad may be contraindicated for certain patient groups, such as lymphoedema patients, and in specific medical situations. In general, heating pads should be avoided in the following cases:

1. **Lymphoedema patients:** Increased swelling: As already mentioned, heat can dilate the blood vessels and worsen the swelling in the affected tissue. The application of a heat pad can increase the accumulation of fluid, which leads to additional discomfort in lymphoedema patients.
2. **Acute inflammation** In acute inflammation, heat can aggravate the inflammatory process. Inflamed tissue is often already swollen, and additional heat can increase this swelling and slow down the healing process.
3. **Fresh injuries or bruises:** In the acute phase following an injury or bruise (usually within the first 48 hours), cold rather than heat is indicated to reduce swelling and bruising. Heat in this case can increase blood flow and make the swelling worse.
4. **Circulatory disorders:** In people with restricted circulation (e.g. arterial vascular disease), the heat pad can cause problems by dilating the blood vessels and having an unfavorable effect on blood flow.
5. **Nerve disorders or loss of sensitivity:** Patients with nerve damage (e.g. diabetes) or loss of sensitivity should use heat pads with caution as they may not feel when the skin is overheated. This increases the risk of burns or skin damage.
6. **Fever:** During a fever, the body is already heated, and additional external heat from a heating pad can further increase body temperature, which is undesirable.

In all these cases, it is advisable to consult a doctor before using a heating pad to rule out any possible risks.

PRODUCT SPECIFICATIONS

Product description	Vibrating Pillow
Model	M260
Power supply	<p>DC motor</p> <p>Input voltage: DC 3.7-5V Rated current: 120mA Rated power: 0.4W Rated speed: 7000 ± (15%) rpm</p> <p>Lithium Polymer Batteries</p> <p>Rated voltage: 3.7V Battery capacity: 3000mAh (3-4 hrs continuous use)</p>
Input power	15W
Safety category	Class II
Packaging size	67.5 x 52 x 41 cm
Product size	66 x 39,5 x 16,5 cm
Massage modes	1-5
Heating levels	35-65°C
Materials	<p>Pillow: 90 % Polyester, 10 % Spandex Foam: Polyurethane</p>

INSTRUCTIONS FOR USE

1. Insert the adapter jack into the DC socket of the pillow. The controller will light up in red, blue and green. The pillow also operates with batteries. If the intensity of the vibrations weakens while using batteries as the power source, switch to using the pillow with the adapter – this simultaneously recharges the batteries.
2. To turn on the heating function, press the **HEATING button**  once. The light will turn on when the power is on. Press and hold the button to turn the heating function off again.
3. To start the oscillating massage in each zone, press the **MASSAGE button** . The light will turn on when massage mode is activated. Press and hold the button to turn the massage function off again.
4. To change modes, press the **MASSAGE button**  until the desired mode is selected. The light will change depending on the selected mode:

Motor	Mode 1			
1	5 seconds			
2	5 seconds			
3		5 seconds		
4		5 seconds		
5			5 seconds	
6			5 seconds	
7				5 seconds
8				5 seconds
Pause	1 second			

Motor	Mode 2			
1	3 seconds	3 seconds	3 seconds	3 seconds
2	3 seconds	3 seconds	3 seconds	3 seconds
3		3 seconds	3 seconds	3 seconds
4		3 seconds	3 seconds	3 seconds
5			3 seconds	3 seconds
6			3 seconds	3 seconds
7				3 seconds
8				3 seconds
Pause	1 second			

Motor	Mode 3			
1	5 seconds		5 seconds	
2	5 seconds		5 seconds	
3	5 seconds		5 seconds	
4	5 seconds		5 seconds	
5		5 seconds		5 seconds
6		5 seconds		5 seconds
7		5 seconds		5 seconds
8		5 seconds		5 seconds
Pause	1 second			

Motor	Mode 4			
1	5 seconds		5 seconds	
2	5 seconds		5 seconds	
3		5 seconds		5 seconds
4		5 seconds		5 seconds
5	5 seconds		5 seconds	
6	5 seconds		5 seconds	
7		5 seconds		5 seconds
8		5 seconds		5 seconds
Pause	1 second			

Motor	Mode 5			
1	5 seconds	5 seconds	5 seconds	5 seconds
2	5 seconds	5 seconds	5 seconds	5 seconds
3	5 seconds	5 seconds	5 seconds	5 seconds
4	5 seconds	5 seconds	5 seconds	5 seconds
5	5 seconds	5 seconds	5 seconds	5 seconds
6	5 seconds	5 seconds	5 seconds	5 seconds
7	5 seconds	5 seconds	5 seconds	5 seconds
8	5 seconds	5 seconds	5 seconds	5 seconds
Pause	1 second			

5. To change the heating level, press the **HEATING button** . The light will switch depending on the selected heating mode:

Heating mode		
Level	Temperatur	Color
1	45°C	Green
2	55°C	Blue
3	65°C	Red


NOTE: Discontinue use immediately if you feel any discomfort or pain.

NOTE: For your benefit, do not use the vibration pillow for more than 60 minutes at a time.

BLUETOOTH FUNCTION

The pillow can be controlled wirelessly via the Venen Engel App: In your app store, search for Venen Engel and download the app to use the Bluetooth feature.

Click on the Venen Engel App, and the device will automatically scan for a connection, the UI instruction are below:

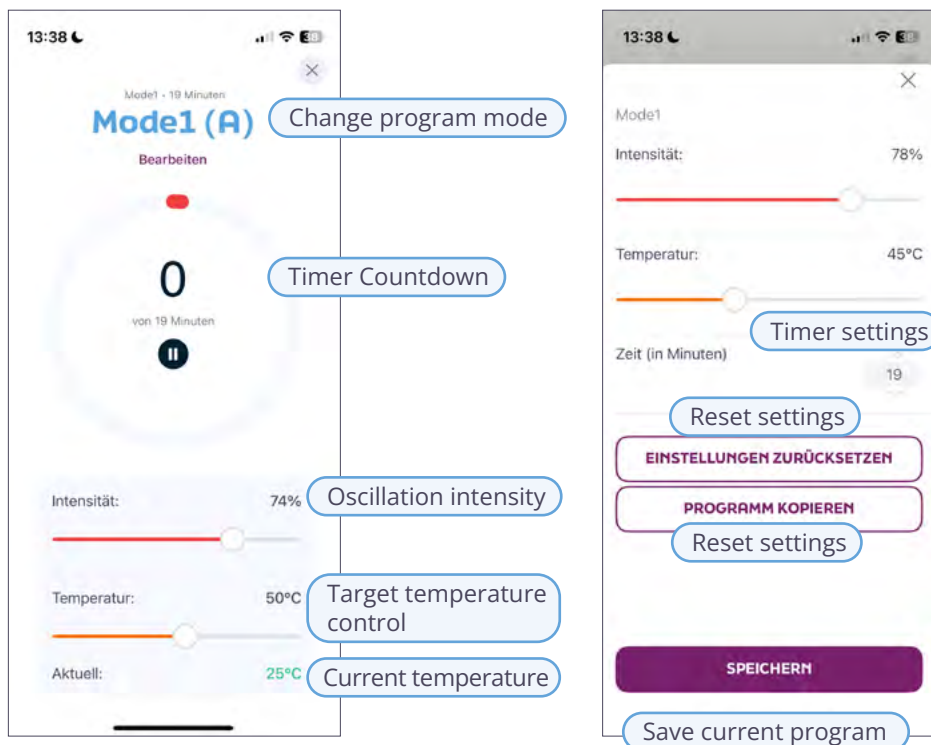Using the Venen Engel app, you can easily change settings such as temperature, duration and intensity of oscillation as well as define and change your own sequences.


In addition, a timer function shows you how long the respective sequence will last. You can also control and customize the duration of the treatment via the app.

USING THE VENEN ENGEL APP


1. Download the Venen Engel App from the AppStore or Google Play Store.
2. Start the app and create an account (if you want to save settings longterm) or continue without an account (in this case, however, your settings will not be saved once you delete the app from your smartphone).
3. Immediately after logging in, the app will attempt to connect to your Vital Flow massage pillow. To do this, allow the Bluetooth function and make sure that your massage pillow is switched on.
4. The My area landing page, displays all the pre-programmed programs (A–E) which you can individually adapt to your preferences when clicking either the program itself or the three dots . If you click on “edit”, you will be able to change the intensity of the oscillation, the temperature, as well as the timer.

You can also reset any changes to the original settings or copy the program to create your own.




To start any of the programs, click on the one you want and then “start”.

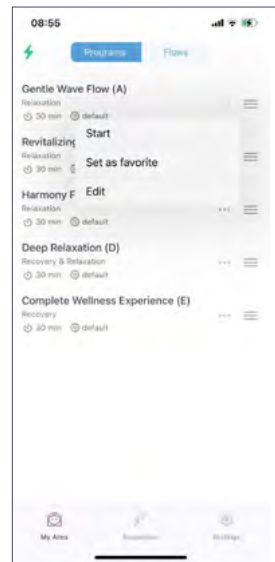5. Also in the **My area** page, you can view the Flows you created (or create your very first own Flow). Just click on “Create Program Flow”, select your preferred program(s) - or change their sequence - and click on “save”.
6. The **Sequences** page provides pre-programed sequences to match specific needs (relaxation, recovery or both) and they, too, can be modified by clicking on the sequence (e.g. “The Gentle Morning Wave”), then the individual Flow (e.g. “Gentle Wave Flow (A)”) and then “Edit Program”.
7. The **Settings** page allows you to change your profile and account settings as well as select or change the device that was previously connected to your phone. The “Help Center” directs you to the Venen Engel website and the Frequently Asked Questions. You can also contact our customer service from there.

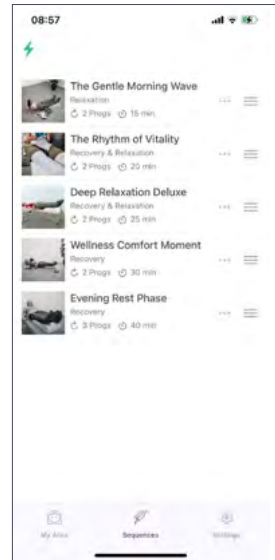HOW TO DELETE YOUR ACCOUNT IN THE VENEN ENGEL APP

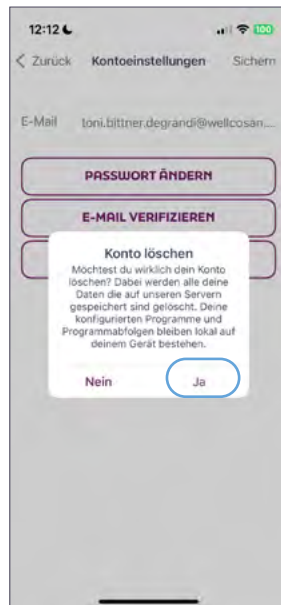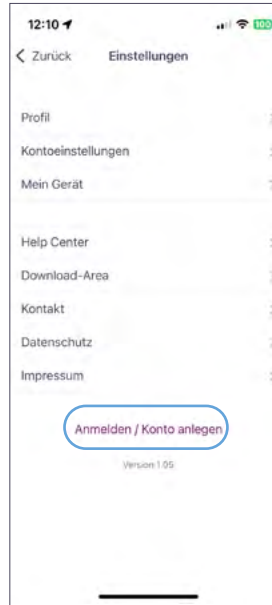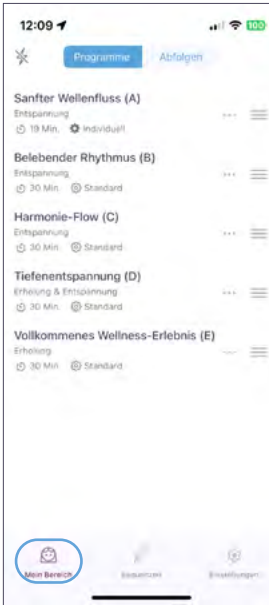MAINTENANCE INSTRUCTIONS

STORAGE:

Place the pillow in its box in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric. Do not wrap the power cord around the pillow. Do not hang the pillow by the power cord.

CLEANING:

- Ensure the pillow is unplugged from the power source and all buttons are off. Allow to cool before cleaning, then wipe down with a damp, clean cloth. Do not use any strong chemical cleaners
- Do not place the pillow in any liquid! Clean the pillow case by hand washing or washing machine after stripping it from the pillow.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass, furniture polish, paint thinner or water hotter than 60°C to clean the pillow case.

SYMBOLS

Symbols	Meaning
	Manufacturer
	Serial number
	Manufacturing date
	Conformity symbol
	Manufacturing country: China
	Product is compliant with 2011/65/EU
	Protection class II (protective insulation)
	Read manual before use
	Not suitable for children under 3 years of age
	No sharp objects
	Do not fold or crumple
	Do not wash
	Do not bleach
	Do not tumble dry
	Hand wash only
	WEEE (electric and electronic used appliance); old appliances must be disposed of according to legal regulations.
	Use only indoors
	Keep dry/protect from water
	France: Disposal instructions for electronic products
	France: Disposal instructions packaging

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Wellcosan GmbH